

TD-1107

EAR INFRARED THERMOMETER




Apply Probe Filter



Measuring



Eject Probe Filter

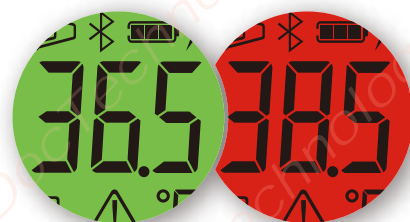

Illuminated Display



iOS

Android



Healthy Check



iOS

Android



Healthy Check Pro

TD-1107

EAR INFRARED THERMOMETER




Probe filter

Hygienic and disposable probe filter for medical applications. Provided with one-hand easy operation of probe filter attachment and ejection without touching probe filter directly to reduce cross contamination.



Illuminated Display

Illuminated display with dual-colored indicator for Hi/Low and fever guidance.

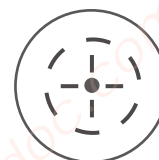

Accurate and Fast

Guiding indicator to ensure proper and stable probe-tip placement to obtain clinical accuracy in seconds. The advanced measuring technology captures and analyzes multiple measurements and selects the most accurate results.



THERMOMETER SPECIFICATIONS

Measurement Range	32.0°C (89.6°F) ~ 43.0°C (109.4°F)
Accuracy	±0.2°C (±0.4°F) for the range of 35.0°C to 42.0°C (95°F to 107.6°F) ±0.3°C (±0.5°F) for the range of < 35.0°C (95°F) or > 42.0°C (107.6°F)
Operating Condition	10°C to 40°C (50°F to 104°F), 30% to 85% RH (non-condensing)
Storage/Transportation Condition	-20°C to 60°C (-4°F to 140°F), 30% to 85% RH (non-condensing)
Resolution	0.1°C / 0.1°F
Measurement Units	°C or °F Switchable
LCD Back Light	<div> <div>32°C</div> <div>38°C</div> <div>43°C</div> </div> <div> <div>Lo</div> <div>Normal</div> <div>Fever</div> <div>Hi</div> </div>
Power source	AA x 2 alkaline battery
Battery life	At Least 3000 Measurements
Memory Capacity	8 Measurements
Power Saving	N/A
Communication	Bluetooth 4.0 (optional)
Dimension	149.56 (L) x 41.7 (W) x 56.15 (H) mm
Weight	80.4g

SHIPPING INFO

Units per Carton	30
Gross Weight per Carton	10.8 kg
Carton Size	40.3 (L) x 33.8 (W) x 40 (H) cm



www.taidoc.com

TaiDoc Technology Corp.


Copyright TaiDoc Technology Corporation. All rights reserved.

TD-1107

EAR INFRARED THERMOMETER

PARTS AND ACCESSORIES



Part No.	Description	Specification	Picture
5507A-001	16-piece probe cover tray	Dimension: 162.34 (L) x 52.84 (W) x 27.4(H) mm	
5507A-002	16-piece probe cover cassette	320 pieces per box 6 boxes per carton Dimension: 148.76 (L) x 46 (W) x 21.7 (H) mm Carton size: 31.5 (L) x 30 (W) x 15 (H) cm Gross weight per carton: 3 KGS	
5507B-001	Cradle	Dimension: 162.7 (L) x 46.7 (W) x 45 (H) mm	
5507B-002	10-piece probe cover cassette	300 pieces per box 30 boxes per carton Dimension: 125 (L) x 86 (W) x 126.5 (H) mm Carton size: 473(L) x 405(W) x 306(H) mm Gross weight per carton: 11.35 KGS	



www.taidoc.com

TaiDoc Technology Corp.

Copyright TaiDoc Technology Corporation. All rights reserved.