


# TD-5086 AKTIVSAFE LANCET



**E**asy to use

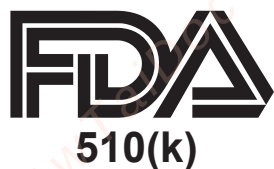

**S**ingle-use only




**F**ive different size




**L**ancet automatic retraction




**F**DA 510(k) cleared Class II medical device

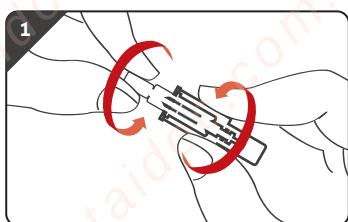

1  
Twist off the protective cap and remove the cap

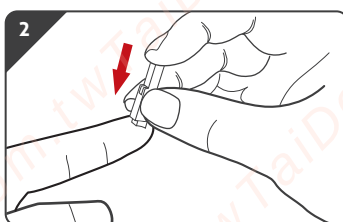

2  
Press the safety lancet body firmly against the puncture site to activate the safety lancet

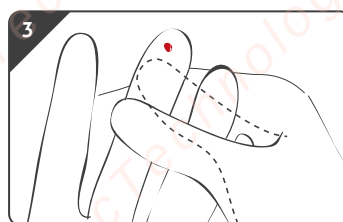

3  
Dispose the used safety lancet and collect blood sample

## GAUGE

21G

23G

26G

28G

30G



|                   |                |               |               |              |              |
|-------------------|----------------|---------------|---------------|--------------|--------------|
| Lancet color      | Green          | Pink          | Orange        | Blue         | Yellow       |
| Penetration Depth | 2.2 mm         | 2.2 mm        | 1.8 mm        | 1.8 mm       | 1.5 mm       |
| Blood Volume      | 60-125 $\mu$ L | 30-75 $\mu$ L | 10-40 $\mu$ L | 5-30 $\mu$ L | 3-20 $\mu$ L |
| Packaging         | 100 pcs/box    | 100 pcs/box   | 100 pcs/box   | 100 pcs/box  | 100 pcs/box  |

## SHIPPING INFO

|                         |                             |
|-------------------------|-----------------------------|
| Units per Carton        | 100                         |
| Gross Weight per Carton | 10.9 kg                     |
| Carton Size             | 90 (L) x 63 (W) x 68 (H) cm |



TaiDoc Technology Corp.

Copyright TaiDoc Technology Corporation. All rights reserved.

www.taidoc.com